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
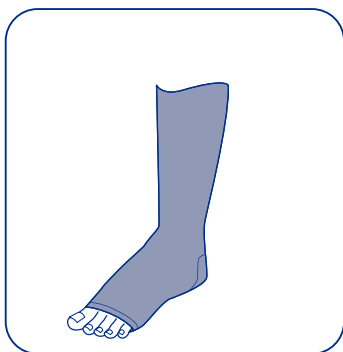
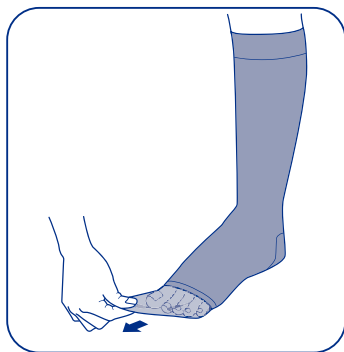
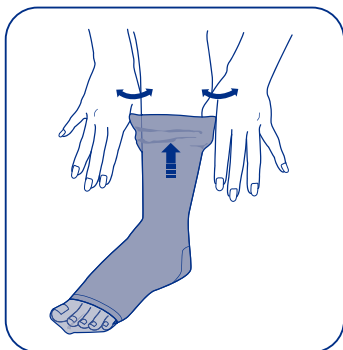
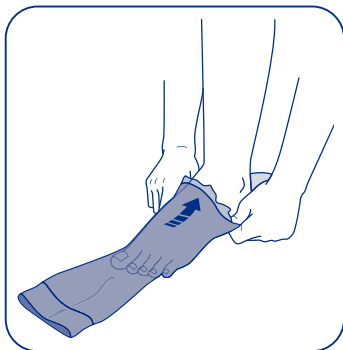
INSTRUCTIONS FOR USE

- IT** Tutore elastico in trama circolare
- EN** Circular weave elastic garment
- FR** Vêtement élastique en trame circulaire
- ES** Prenda elástica de trama circular
- DE** Rundgestrickte elastische Kompressionshilfe
- PT** Tutor elástico de trama circular
- PL** Wyrób uciskowy o splocie kolistym
- DA** Rundvævet elastisk støtteartikel
- NO** Sirkulær vevd elastisk skinne
- RU** Компрессионное трикотажное изделие кругового плетения
- NL** Elastische rondbreibrace
- HR** Elastični steznik kružne potke
- EL** Ελαστικά προϊόντα διαβαθμισμένης συμπίεσης κυκλικής πλέξης
- AR** دعامة مرنة من نسيج دائري
- HI** परिपत्र बुनाई में इलास्टिक ब्रेस
- TR** Dairesel dokuma elastik korse

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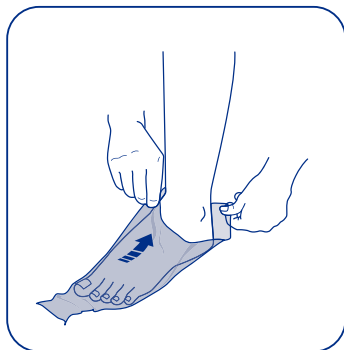Circular weave elastic-compression garments manufactured in Italy in compliance with the main international product and service quality standards. Products periodically tested and certified by international reference bodies.

RAL COMPRESSION CLASSES

- I class** (18-21 mmHg)
- II class** (23-32 mmHg)
- III class** (34-46 mmHg)
- IV class** (>49 mmHg)

ASQUAL COMPRESSION CLASSES

- I class** (10-15 mmHg)
- II class** (15-20 mmHg)
- III class** (20-36 mmHg)
- IV class** (>36 mmHg)

GENERAL INDICATIONS TO PUT ON CIZETA MEDICALI ELASTIC GARMENTS CORRECTLY

1. Compression treatments with a medical elastic garment should be carried out only during daytime.
2. Before putting on the elastic garment, it is advisable to keep the arm or the leg in a higher position than the torso for at least 5 minutes, to facilitate blood reflux. It would be perfect to put on the elastic garment in the morning before getting up.
3. Put on the elastic garment on clean and dry limbs.
4. In case the use of creams and/or ointments for the treatment of the limb is required, it is preferable to apply them in the evening. In the morning, wash the affected limb carefully before putting on the device. Oils, skin ointments, etc. can affect the effectiveness and durability of the device.
5. It is suggested to use rubber gloves to prevent rings or nails from damaging the device.

INSTRUCTIONS TO PUT ON THE ARMBANDS CORRECTLY

1. Turn the armband inside out.
2. Insert your hand and slide the armband along the limb, by accompanying it gently without pulling it on strongly.

INSTRUCTIONS TO PUT ON THE STOCKINGS CORRECTLY

1. Turn the stocking inside out up to the heel.
2. Insert the foot, put the stocking heel in the

correct position and then make the stocking slide along the limb, by accompanying it gently without pulling it on strongly.

3. To put on the stocking, the VARITEC® stocking wearing aid or the VARISLIDE® inserter can be helpful.

INTENDED USE

The intended use of CIZETA MEDICALI elastic-compression devices is specific for the treatment of venous and lymphatic system disorders and in case of lipoedema, namely a disorder of adipose tissue depending on oestrogens.

TARGET PATIENTS

CIZETA MEDICALI elastic-compression devices are used in the treatment of people suffering from phlebopathies and lymphopathies and lipoedema after clinical evaluation by the competent physician.

INDICATIONS

- Functional disorders, such as feeling of heaviness and tiredness of the legs in presence of varicose veins.
- Evening localised oedema (of venous origin).
- Post-operative and post-traumatic oedema.
- Post-thrombotic oedema with or without secondary treatment of varices.
- Deep and Superficial Venous Thrombosis (DVT and SVT).
- Chronic Venous Insufficiency (CVI) in all evolving clinical levels (symptoms of venous insufficiency without signs, reticular varices and telangiectases, varicose veins, venous oedema, trophic skin disorders due to CVI, venous ulcers in their active phase and for the prevention of recurrences).
- Secondary venous stasis syndrome and arthropathies.
- Treatment following:
 - sclerosing therapy
 - selective varicose vein surgery
 - stripping
- Reversible lymphedema.
- Lipoedema.
- Angiodysplasia (in selected cases and with specific braces).
- Venous disease prophylaxis in patients at risk.
- Prophylaxis of thrombosis and embolism.
- Chronic venous stasis ulcers (treatment

and prophylaxis of recurrences).

- Varicose veins during pregnancy without oedema.
- Varicose veins during pregnancy with oedema.
- Primary and secondary varicose veins.

CONTRAINDICATIONS

- Peripheral obliterative arteriopathy with ankle/arm index below 0,7 (ABI) or peripheral blood pressure below 80 mmHg or uncompensated.
- Immobility/severe hypomobility (except for A.T.E.[®] anti-thrombo embolic stockings).
- Acute skin infections.
- Exudative dermatitis.
- Acute dermal hypo-dermatitis.
- Diseases of the lower limbs other than those of venous-lymphatic origin.
- Oedemas not having venous origin (dysmetabolic, diabetic, renal, cardiac, endocrine...) in the acute phase or undiagnosed.
- Known history of contact allergy caused by medical compression stockings and their components.
- Congenital/acquired anatomical anomalies (except for the possibility of using custom-made elastic garments).

PRODUCT STORAGE

Store in a cool and dry place, and in its original packaging. Protect from direct sunlight.

Washing instructions

1. Machine wash at a temperature not exceeding 40°C.
2. Preferably use neutral soap.
3. Run a gentle spin cycle.
4. Dry on flat surfaces away from heat sources, preferably not in the sun.
5. Remember that frequent washing does not damage the stocking, instead it prolongs its life.

MANAGEMENT OF DAMAGED DEVICES

If the device is damaged, notify the manufacturer immediately.

SIDE EFFECTS AND SERIOUS ACCIDENTS

Currently, there are no known side effects or serious accidents. Should negative episodes arise during use, consult your doctor. In the event of a serious accident involving the device, notify the manufacturer and the competent doctor immediately.

COMPOSITION

The composition is indicated inside (label or package).

VALIDITY OF CIZETA MEDICALI COMPRESSION BRACES

The high-quality materials used to manufacture the elastic-compression garments, the technologies adopted in its technical structure and the product certifications achieved suggest a 6-month validity for the effectiveness of the elastic garments, starting from their first use. The need for a possible replacement in periods prior to said deadline may arise, if the instructions indicated in this leaflet are not observed.

REMARKS

The level of compression and the type of medical elastic garment should be suggested by the doctor. The evolution of the therapy or changes of the clinical conditions could result in a reduction of the oedema and therefore require the prescription of a new elastic garment.

	CE mark
	Manufacturer
	Date of production
	Expiry date
	Medical device
	Catalogue number
	Lot code
	Unique Device Identification
	Do not use if package is damaged
	Wash at 40°C
	Extra delicate wash
	Horizontal drying

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CIZETA MEDICALI S.p.A.