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INSTRUCTIONS FOR USE AND COUNTRY-SPECIFIC NOTES

1. Product description

Medical compression stockings

2. Intended purpose

Circular-knitted and flat-knitted medical compression stockings are designed to exert external, physical compression.

Every JOBST® medical compression stocking has a specific compression level designed to help manage specific venous and lymphatic conditions.

The required compression class and style must be determined by your healthcare professional. The compression level is printed on the package label.

Indications

A medical compression stocking is used for the management of the following indications:

Compression Class 1 (CCL 1)

- Tired, heavy, aching legs
- Minor varices without edema
- Thrombosis prophylaxis in mobile patients
- Mild swelling of the feet, ankles and legs
- Onset of pregnancy-related varices
- Controlled cardiac edema
- Arthritic patients with reasonable dexterity

- Primary or secondary lymphedema in International Society of Lymphology (ISL) stages I-II with no or minimal shape distortion (or lymphedema ISL stage III when higher compression is not tolerated)

Compression Class 2 (CCL 2)

- Tired, heavy, aching legs
- Moderate to severe varices with and without edema
- Moderate to severe pregnancy-related varices with and without edema
- Moderate edema
- Chronic venous insufficiency (CVI) C3-C5 according to CEAP: venous edema, skin alterations and healed ulcer cruris
- Helps to prevent the recurrence of venous ulcers
- For adjunct management of superficial venous thrombosis
- Post vein surgery
- Post sclerotherapy
- Post-traumatic and post-operative edema
- Arthritic patients with reasonable dexterity
- Primary or secondary lymphedema in ISL stage II-III with some shape distortions
- Phlebolympheidema
- Lipedema

Compression Class 3 (CCL 3)

- Severe varices with or without edema
- Severe edema
- Chronic venous insufficiency (CVI) C3-C6 according to CEAP: venous edema, skin alterations, ulcer cruris and healed ulcer cruris
- Helps to prevent the recurrence of venous ulcers
- Post-thrombotic syndrome (PTS)
- Post vein surgery
- Post sclerotherapy
- Active patients with lymphedema in ISL stage III
- Risk of edema rebound
- Severe lipedema

Compression Class 4 (CCL 4)

- Severe post-thrombotic syndrome (PTS)
- Severe (complex) lymphoedema, ISL stage III
- Severe edema of non-arterial origin
- Swelling not contained by medium or high pressure stocking

Please contact your healthcare professional if you are unsure whether the product is appropriate for you.

3. Contraindications

A JOBST® medical compression stocking is contraindicated (should NOT be worn) if you have one or more of the following medical conditions:

- Severe peripheral arterial occlusive disease
- Uncontrolled / decompensated heart failure
- Septic phlebitis
- Phlegmasia coerulea dolens
- Untreated and / or advancing skin infection (e.g. erysipelas) involving the area where the garment will be applied
- Open wounds not covered with an appropriate dressing
- Incompatibility to material

4. Warnings and precautions

A medical compression stocking should only be used after consultation with your healthcare professional and checked regularly by him/her on the appropriateness of your therapy.

Your healthcare professional decides whether compression in general can be applied to patients suffering from polyarthritides (unspecified), complex regional pain syndrome (CRPS, Morbus Sudeck), malignant lymphedema, and medical conditions in which increased venous and lymphatic return is not desired.

Consult your healthcare professional if you have experienced or develop: controlled / compensated heart failure, mild to moderate peripheral arterial occlusive disease, peripheral neuropathy, severely compromised skin sensibility and impaired sensitivity of the limb, as well as skin infection or weeping dermatoses.

If you notice discomfort, pain, numbness, changes in skin color, constrictions, pressure marks, appearance of new wounds or swelling while wearing a compression stocking, stop using the product and contact your healthcare professional. Patient measurements should be taken by a trained person to determine correct product fit.

Incorrectly measured and badly fitting compression products may cause constriction leading to e.g. skin necrosis, pressure ulcers, and pressure-related peripheral nerve damage during use.

When re-ordering a product we recommend remeasuring the patient as the limb circumference might have changed.

Appropriate hygiene and skin care should be applied.

Appropriate footwear must be worn when walking to help prevent falls.

Not recommended to wear during bed rest.

Caution is required if you are confined to bed or unable to walk.

5. Directions for use

To get the best medical result from your JOBST® medical compression stocking, wear it regularly and consistently.

INSTRUCTIONS FOR USE AND COUNTRY-SPECIFIC NOTES

EN

Ideally you should put on your compression stocking in the morning, immediately after getting up. You may experience some difficulty putting on the stocking in the beginning, but this will become easier after a number of days.

When applying the compression stocking, please do not pull on the edge of the stocking instead follow these recommendations below.

Turn the compression stocking inside out all the way to the heel.

Pull apart the folded-in foot part using both hands and pull over the foot all the way to the heel.

Gradually pull the compression stocking over the instep and heel.

Then gradually pull the compression stocking upward by returning the stocking to a right side out position piece-by-piece.

Then distribute the stocking evenly by massaging with both hands until the stocking fits and does not form creases. It is particularly helpful to wear rubber gloves for this step.

For a waist high garment continue working the garment up to the middle of your thigh. Stand up. With both hands, grab a handful of fabric and pull it up to your waist.

Adjust your support so that it is snug at your midsection. This will keep it from sliding down. Pull the midsection to its proper height.

Slip your hands, palm out, into the back and lift the fabric up and away from your skin. Let it fall back to fit snugly over the buttocks.

Special notes for donning/doffing:

Long fingernails and jewelry may damage the stocking when putting them on.

For easier application use rubber gloves or a donning aid.

Do not bring the compression stocking into contact with greasy ointments, body lotions or similar substances unless recommended by JOBST®.

Loose threads or loops should not be cut off, since this may cause holes or runs to form in the stocking.

6. Cleansing and care

We recommend to wash the JOBST® compression stocking on a daily basis, or at least every other day.

For this purpose, you may use a mild detergent and a low spin cycle setting. Do not use laundry additives such as fabric softener, optical brightener or stain remover.

Rinse your compression stockings well.

Roll up the wet compression stockings in a towel after washing and press out excess moisture. Do not wring.

Do not dry the compression stockings on a radiator or in direct sunlight.

Machine washable in warm water.


Do not bleach.

Flat drying or tumble dry – low heat. For tumble dryer compatibility please check your compression stocking's product label; the label should not be removed.

Do not dry clean.

Do not iron.

Note: a laundry net provides additional protection for the knitted fabric.

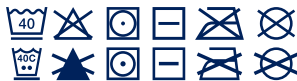Silicone band care:

When needed, hand wash silicone band to remove excess creams, lotions, hair and other foreign material.

7. Storage

A maximum storage period of up to 36 months applies for ready-to-wear products. This is followed by a maximum period of use of 6 months for the compression stocking. The use-by date of the medical device is printed on the package label with the symbol of an hourglass.

Made-to-measure products are products intended for immediate use and are marked with a manufacturing date.

Due to wear or slackening of the material, the stocking is intended for a wearing period of six months. This assumes correct handling (e.g. product care, donning and doffing).

Please protect your compression stocking from direct exposure to sunlight, heat and moisture. We recommend storage at room temperature.

8. Notice

Any serious incident that has occurred in relation to this device should be reported to BSN-JOBST® GmbH and the competent authority of your state. Caution applicable for made-to-measure products: Federal Law (U.S.A.) restricts this device to sale by or on the order of a physician.

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	Artwork ID: JBR14543.05	Review Loop: 02	Printing Colours:  PANTONE 281 C
Project Title: C2847_JOBST_MLW_EMM		Date of Review Loop: 26.06.2023 GT_KLND	
Contact ArtworkManagementMedical@essity.com		Replaces Artwork ID: JBR14543.04	
Contact Address: BSN medical GmbH Schützenstraße 1–3 22761 Hamburg Germany		Technical Drawing for print ID	Guiding Colours:  Diecut - DO NOT PRINT
Languages: EN, DE, FR, ES, PT, IT, NL, SV, NO, DA, FI, IS, ET, LT, LV, PL, HU, RO, CS, SL, SK, RU, TR, AR, HE			

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BSN-JOBST GmbH
Essity UK Ltd, Southfields Road, DUNSTABLE, Beds, LU6 3EJ
BSN medical (Aust.) Pty Ltd